



FRAXEL FTX[®] SKIN RESURFACING, REIMAGINED.

A DEMONSTRATED SAFETY PROFILE

20+

Years of
innovation

50+

Clinical
studies

9+

Active
patents*

1.5+

Million treatments
performed

*2026 data to date on total number of patents.

MODEL

TWO WAVELENGTHS. ONE INTELLIGENT APPROACH.

Fraxel FTX® is a highly versatile and recognised treatment in aesthetics, offering powerful results. It stands out for its precision, yet its true impact lies in its ability to address 5 common skin concerns contributing to skin health

1. Fine lines and wrinkles
2. Lentigos & ephelides
3. Acne scars & surgical scars
4. Dyschromia, including melasma
5. Actinic keratosis

2026 NewBeauty Reader's Choice
favourite laser for sun damage award.



ONE OF A KIND TECHNOLOGY

HOW FRAXEL® WORKS


Consistency and precision with Intelligent Optical Tracking™ and AccuTRAC™.




The Fraxel FTX has two wavelengths to enable targeted treatment of different tissue depths. The 1927 nanometer wavelength targets the epidermal layer. The 1550 nanometer wavelength is able to penetrate into the deeper dermis, while still treating the epidermis.



As the tip moves over the skin, the laser treats a fraction of the skin at a time creating a large number of microscopic treatment zones on the skin in a dot-like pattern.

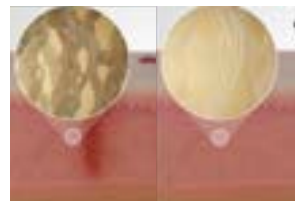

Since the laser leaves the surrounding tissue intact, it promotes rapid healing. It works by targeting the water that's naturally occurring in the skin, stimulating the body's own natural healing process, promoting collagen remodelling and replacing the old and damaged cells with fresh, glowing, healthy-looking skin.

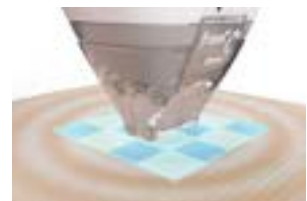

The patented Fast Tracking Xperience with AccuTRAC® optical system calculates hand motion during treatment to help ensure proper contact with skin, dosimetry and uniform delivery. It also makes laser firing more efficient.

FEATURES AND BENEFITS

Console



A modern design with an integrated arm and cable assembly, that is in line with the Solta product family.

- **95%** of providers noted the Fraxel FTX system was easy to use*.

Handpiece (HP)

- Ergonomic HP with **20%** circumference & weight reduction, designed for improved field of view via tapered tip, and treatment precision
- Direct airflow to the skin designed for patient comfort
 - Cooling is routed within HP to the tip
 - Results in **11%** decrease in skin temperature compared to Fraxel DUAL**
 - Airflow is directed at treatment area
- **96%** of providers noted that ergonomic adjustments to reduce fatigue were an improvement to the DUAL system*



Single Use Consumable Tips (15mm & 7mm)



- No need to change rollers, designed to save time
- Reduced risk of patient cross contamination
- Designed for patient comfort, the integrated cooling flows through the handpiece and tip directly to the skin
- Quickly switch between tips to provide the ability to create a full face "360-degree" treatment
- **95%** of providers thought that single-use tips likely decreased the chance of cross-contamination*

Credit Key



- Enhanced security features anti-counterfeit software
- Key provides lockout feature for safety
- Pre-loaded with up to 12 full face treatments
- Credit key allows changing 15mm & 7mm tips within a single treatment
- **77%** of providers reported switching between the 15 and 7 mm tips was faster with the Fraxel FTX compared to Fraxel DUAL*

Refreshed Treatment Screen



Redesigned guided user Interface (GUI) designed to facilitate treatments including:

- Treatment Settings
- Total & estimated energy
- Handpiece Speed
- Colour change between wavelengths

*HCP feedback based on Fraxel FTX Beta Site HCP Survey Results 2025, N=22

** Based on internal preclinical (swine model) testing, demonstrates an 11% decrease in skin temperature compared to Fraxel DUAL. Internal preclinical R&D data on file.



FRAXEL FTX® ALL SKIN TYPES ALL OVER THE BODY

Latest advanced resurfacing technology from Solta.
The Fraxel FTX builds on a legacy of award-winning performance designed to enhance speed and ease-of-use.

- Fraxel® lasers can treat women and men of all ages and skin colours
- Fraxel® can be used on the face and body and works on delicate skin areas like the neck, chest and hands

POWERFUL, FAST-TRACKING MOVEMENT

- 1550 nm wavelength is an infrared Erbium fibre laser that can penetrate into the deeper dermis
- 1927 nm wavelength is an infrared Thulium fibre laser has a shallower depth of penetration, targeting the epidermal layer
- Intelligent optical tracing system with AccuTRAC® provides faster and more accurate delivery of microscopic treatment zones (MTZ). Half of providers surveyed noted the fast tracking system was up to **20% faster** than Fraxel DUAL.*

*HCP feedback based on Fraxel FTX Beta Site HCP Survey Results 2025, N=22

A PROFITABLE TREATMENT FOR YOUR PRACTICE*

- Patients ask for Fraxel® laser treatment by name
- Multiple treatment opportunities per patient
- Treatment can be delegated in certain jurisdictions**
- Fraxel FTX has a demonstrated safety profile but may not be right for everyone.

*Solta does not guarantee business success.

**Check your local jurisdictions for delegable.



DYSCHROMIA INCL. MELASMA

- Targets ageing and sun-damaged skin to help improve tone, texture, and radiance
- Helps improve the clinical appearance of photodamage and photoageing



Before



Before

Photographs courtesy of Dr Brian Biesman, MD



After

2 Weeks After 3 Treatments*



After

4 Weeks After 1 Treatment*

ACNE SCARS & SURGICAL SCARS

- Helps improve the appearance of acne scars and surgical scars



Before

1 Day After 1 Treatment*



Before

Photographs courtesy of Dr. Ophelia Veraitch



After

1.5 Months After 4 Treatments*



After

3 Months After 5 Treatments*

*Number of treatments is at the discretion of the treating provider.
Unretouched photos. Individual results may vary.

WRINKLES & FINE LINES

- Fraxel® helps rejuvenate skin to help smooth wrinkles and fine lines. Remodels deep dermal collagen and turns over epidermal tissue, while sparing the stratum corneum.



Before

Photographs courtesy of Dr Suzanne L Kilmer, M.D.



Before

Photographs courtesy of Steven Weiner, M.D.



After

4 Weeks After 1 Treatment*



After

1 Month After 3 Treatments*

ACTINIC KERATOSIS



Before

Photographs courtesy of Dr Jill Waibel, MD



Before

Photographs courtesy of Dr Jill Waibel, MD



After

4 Weeks After 1 Treatment*



After

4 Weeks After 1 Treatment*

*Number of treatments is at the discretion of the treating provider.
Unretouched photos. Individual results may vary.

This product is not available for purchase by the general public. Always follow the directions for use.
FRAXEL® is a trademark of Solta Medical, Inc. or its affiliates ©2026 Solta Medical, Inc. or its affiliates.
FRXFTX.0015.AU.26

SOLTA SUPPORT TO DRIVE SUCCESS

TRAINING SUPPORT

Solta Academy is an online educational platform that provides opportunities to physicians and their staff to build on their clinical knowledge and gain insights at their own pace and time to help deliver exceptional patient care.

This is followed by practical application training with our Clinical Training Specialists.

PRACTICE SUPPORT

Offering assistance with business planning, marketing strategies and support. Includes placement on Solta Medical's clinic finder site, customised bespoke marketing materials alongside a full suite of ongoing complimentary digital and print resources, ensuring your marketing efforts print resources are always effective.

PRODUCT SUPPORT

Standard 2 -year manufacturer warranties with access to responsive Product Support Team that includes field service engineers for assistance in minimising business disruption.



SOLTA MEDICAL®

FOLLOW US ON

Solta Medical Australia

Fraxelau

Fraxel AU

www.Fraxel.com.au